

gingifast



Technical
A-Silicone for Gingival Mask

Zhermack 
Dental

gingifast elastic

Natural aesthetics.

APPLICATIONS

Fast gingival reproduction in A-Silicone delivered in a **self-mixing** system. Allows the correct morphology in the design of fixed prostheses.

AESTHETICS

The **aesthetic result** of **GINGIFAST ELASTIC** is obtained from its **translucency** and veining effect which give it a natural appearance.



1 Master model with a mask made in Zetalabor Lab Putty



2 Application of Gingifast Separator on to the putty mask obtained with Zetalabor



3 Space ground away on the master model ready for the injection of Gingifast Elastic



4 Positioning of the mask in Zetalabor on the model and injection of Gingifast Elastic



5 Detail of the reproduced tissue after finishing



6 Work finished on the master model

gingifast

gingifast rigid


Versatile and easy to trim.

APPLICATIONS

The gingival reproduction used in **implant prosthesis simulates** the conditioning of the gingival tissue for the determination of the shape of the cervical contours and the emergence profile of the prosthesis.

HARDNESS


The **balanced hardness** makes **trimming** of the edges easier when finishing the gingival tissue on the model.

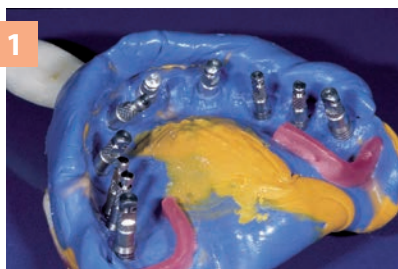


1 Boxing the impression with wax



2 Application of Gingifast Separator onto the impression

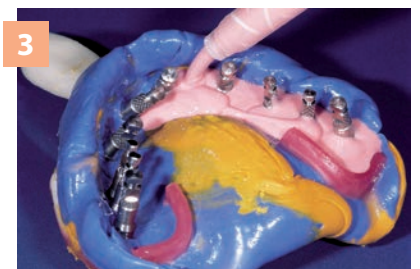

3 Immersed mixing tip during the extrusion for extra detail and homogenous distribution

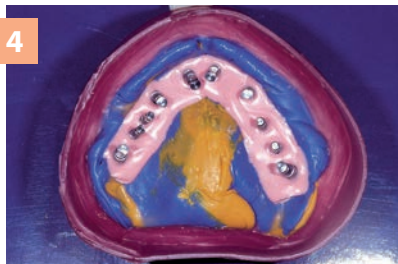

4 Boxed impression ready for final casting of the model



5 Detail of the reproduced tissue



6 Easy trimming

gingifast



Packaging

C401500 GINGIFAST ELASTIC	2 x 50 ml cartridges +12 x oral tips, yellow + 1 x Separator 10 ml + 12 x oral tips, yellow + cap spray
C401520 GINGIFAST RIGID	2 x 50 ml cartridges +12 x oral tips, yellow + 1 x Separator 10 ml + 12 x oral tips, yellow + cap spray

Accessories

C202100 DISPENSER D2	Dispenser D2 - 1:1
C202070 MIXING TIPS SMALL	Mixing tips, yellow (48 pcs)
C202090 ORAL TIPS	Intraoral tips (48 pcs)
C400888 SEPARATOR FOR GINGIFAST	2 bottles Separator, 10 ml

BRANCHES

GERMANY

Zhermack GmbH Deutschland
Öhlmühle 10 - D-49448 Marl

T +49 (0) 54 43 / 20 33 - 0
F +49 (0) 54 43 / 20 33 - 11
info@zhermack.de
www.zhermack.de

POLAND

Zhermapol Sp. z o.o.
ul. Augustówka 14
02 - 981 Warszawa - Polska

T +48 22 858 82 72
+48 22 292 71 17
F +48 22 642 07 14
biuro@zhermack.com
www.zhermapol.pl

USA

Zhermack Inc.
P.O. Box 4195, River Edge
New Jersey 07661 - 4195

T +1 (877) 819 6206
+1 (732) 389 8540
F +1 (732) 389 8543
info@zhermackusa.com
www.zhermackusa.com