

ELITE MODEL
THIXOTROPIC

GESSO TIPO 3 PER MODELLI IN GENERE E ANTAGONISTI
TYPE 3 DENTAL STONE FOR ANTAGONISTS AND MODELS IN GENERAL
HARTGIPS, KLASSE 3, FÜR SITUATIONS- UND GEGENBISSMODELLE
PLÂTRE TYPE 3 POUR MODÈLES EN GÉNÉRAL ET D'ANTAGONISTES
YESO TIPO 3 PARA MODELOS EN GENERAL Y ANTAGONISTAS


WATER / POWDER RATIO	30 ml / 100 g
WORKING TIME	7'
SETTING TIME (VICAT)*	12'
SETTING EXPANSION 2h*	0.04 %
COMPRESSIVE STRENGTH 1h*	31 MPa

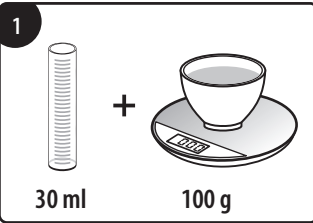
Rapporto acqua/polvere
Water/powder ratio
Wasser/Pulver Verhältnis
Ratio eau/poudre
Dosificación agua/polvo

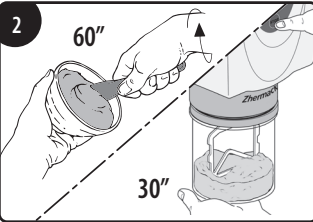Miscelare per 60/30 secondi
Mix for 60/30 seconds
60/30 Sekunden anmischen
Mélanger pendant 60/30 secondes
Mezclar entro 60/30 segundos

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


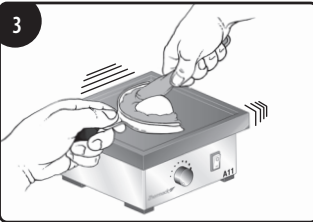
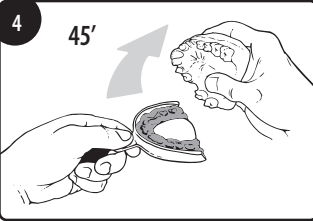
WATER / POWDER RATIO	30 ml / 100 g
WORKING TIME	4'
SETTING TIME (VICAT)*	8'
SETTING EXPANSION 2h*	0.1 %
COMPRESSIVE STRENGTH 1h*	36 MPa

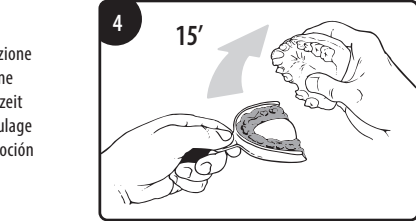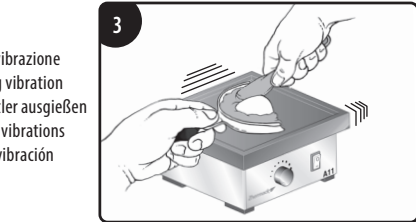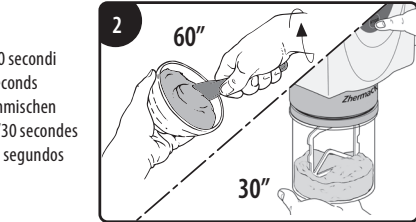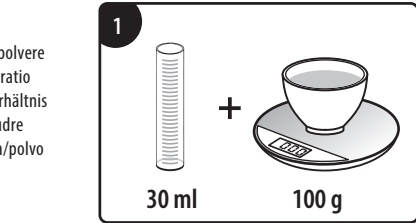
ELITE ARTI
THIXOTROPIC

GESSO TIPO 3 PER ARTICOLATORI
TYPE 3 ARTICULATING DENTAL STONE
ARTIKULATIONSGIPS, KLASSE 3
PLÂTRE TYPE 3 POUR ARTICULATEURS
YESO TIPO 3 PARA ARTICULADORES


WATER / POWDER RATIO	30 ml / 100 g
WORKING TIME	3'
SETTING TIME (VICAT)*	5'
SETTING EXPANSION 2h*	0.02 %
COMPRESSIVE STRENGTH 1h*	26 MPa

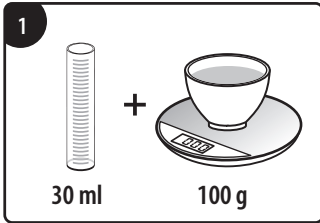Rapporto acqua/polvere
Water/powder ratio
Wasser/Pulver Verhältnis
Ratio eau/poudre
Dosificación agua/polvo

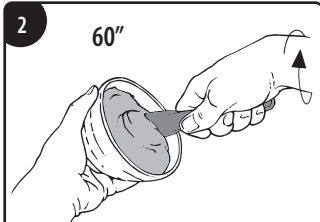Miscelare per 60 secondi
Mix for 60 seconds
60 Sekunden anmischen
Mélanger pendant 60 secondes
Mezclar entro 60 segundos


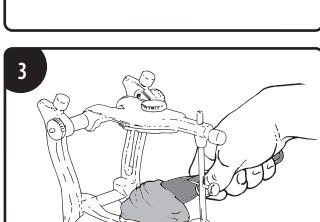
Applicare sull'articolatore
Put on articulator
In den Artikulator einsetzen
Appliquer sur l'articulateur
Aplicar en el articulador


ELITE ARTI FAST
THIXOTROPIC

GESSO RAPIDO TIPO 3 PER ARTICOLATORI
TYPE 3 FAST SETTING ARTICULATING DENTAL STONE
SCHNELLABBINDENDER ARTIKULATIONSGIPS, KLASSE 3
PLÂTRE RAPIDE TYPE 3 POUR ARTICULATEURS
YESO RÁPIDO TIPO 3 PARA ARTICULADORES

WATER / POWDER RATIO	30 ml / 100 g
WORKING TIME	2'
SETTING TIME (VICAT)*	3'
SETTING EXPANSION 2h*	0.02 %
COMPRESSIVE STRENGTH 1h*	26 MPa

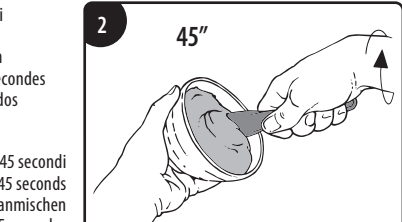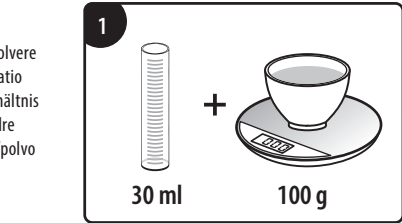
Miscelare per 45 secondi
Mix for 45 seconds
45 Sekunden anmischen
Mélanger pendant 45 secondes
Mezclar entro 45 segundos

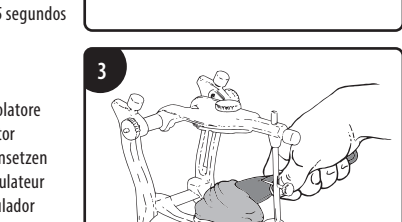
ELITE ORTHO
THIXOTROPIC

GESSO TIPO 3 PER MODELLI IN ORTODONZIA
TYPE 3 DENTAL STONE FOR ORTHODONTIC MODELS
KFO GIPS, KLASSE 3
PLÂTRE TYPE 3 POUR MODÈLES EN ORTHODONTIE
YESO TIPO 3 PARA MODELOS, EN ORTODONCIA


WATER / POWDER RATIO	24 ml / 100 g
WORKING TIME	7'
SETTING TIME (VICAT)*	14'
SETTING EXPANSION 2h*	0.04 %
COMPRESSIVE STRENGTH 1h*	32 MPa

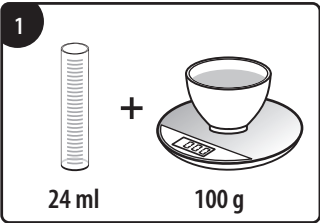
Rapporto acqua/polvere
Water/powder ratio
Wasser/Pulver Verhältnis
Ratio eau/poudre
Dosificación agua/polvo

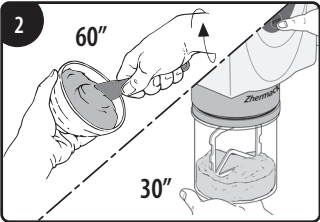Miscelare per 60/30 secondi
Mix for 60/30 seconds
60/30 Sekunden anmischen
Mélanger pendant 60/30 secondes
Mezclar entro 60/30 segundos

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
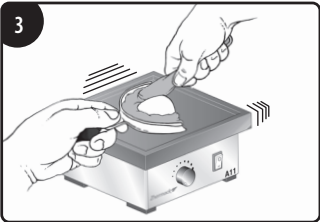
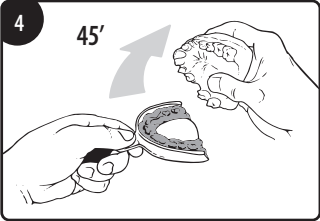
WATER / POWDER RATIO	21 ml / 100 g
WORKING TIME	12'
SETTING TIME (VICAT)*	14'
SETTING EXPANSION 2h*	0.08 %
COMPRESSIVE STRENGTH 1h*	50 MPa

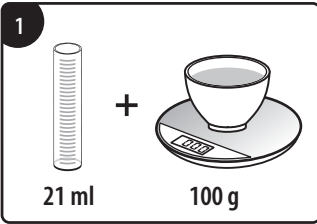
Rapporto acqua/polvere
Water/powder ratio
Wasser/Pulver Verhältnis
Ratio eau/poudre
Dosificación agua/polvo

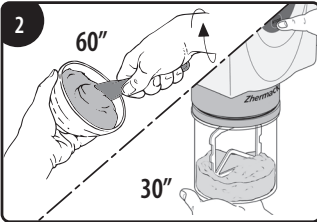Miscelare per 60/30 secondi
Mix for 60/30 seconds
60/30 Sekunden anmischen
Mélanger pendant 60/30 secondes
Mezclar entro 60/30 segundos

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
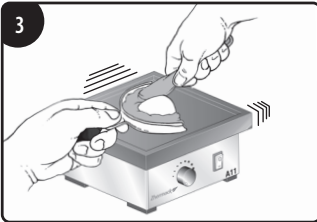
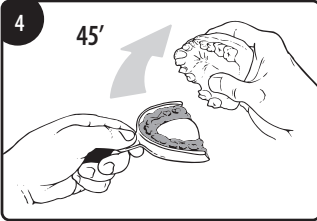
WATER / POWDER RATIO	20 ml / 100 g
WORKING TIME	12'
SETTING TIME (VICAT)*	14'
SETTING EXPANSION 2h*	0.08 %
COMPRESSIVE STRENGTH 1h*	52 MPa

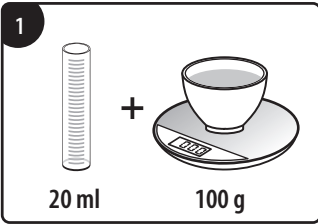
Rapporto acqua/polvere
Water/powder ratio
Wasser/Pulver Verhältnis
Ratio eau/poudre
Dosificación agua/polvo

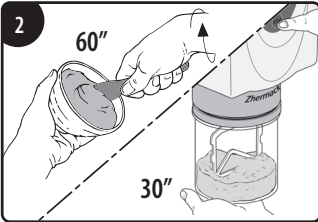Miscelare per 60/30 secondi
Mix for 60/30 seconds
60/30 Sekunden anmischen
Mélanger pendant 60/30 secondes
Mezclar entro 60/30 segundos

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
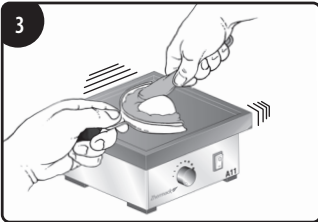
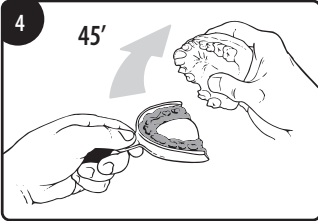
WATER / POWDER RATIO	20 ml / 100 g
WORKING TIME	5'
SETTING TIME (VICAT)*	9'
SETTING EXPANSION 2h*	0.08 %
COMPRESSIVE STRENGTH 1h*	52 MPa

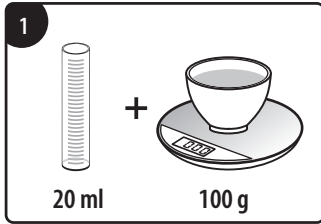
Rapporto acqua/polvere
Water/powder ratio
Wasser/Pulver Verhältnis
Ratio eau/poudre
Dosificación agua/polvo

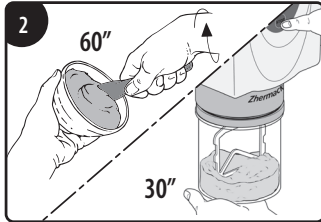Miscelare per 60/30 secondi
Mix for 60/30 seconds
60/30 Sekunden anmischen
Mélanger pendant 60/30 secondes
Mezclar entro 60/30 segundos

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
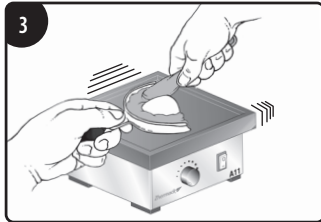
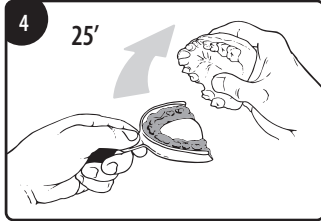
WATER / POWDER RATIO	25 ml / 100 g
WORKING TIME	8'
SETTING TIME (VICAT)*	14'
SETTING EXPANSION 2h*	0.08 %
COMPRESSIVE STRENGTH 1h*	42 MPa

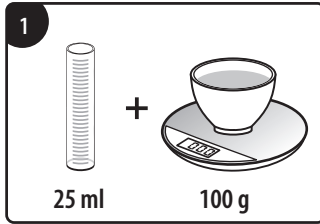
Rapporto acqua/polvere
Water/powder ratio
Wasser/Pulver Verhältnis
Ratio eau/poudre
Dosificación agua/polvo

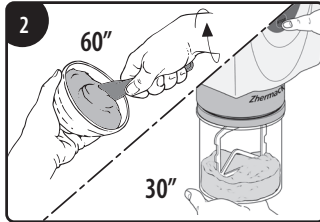Miscelare per 60/30 secondi
Mix for 60/30 seconds
60/30 Sekunden anmischen
Mélanger pendant 60/30 secondes
Mezclar entro 60/30 segundos

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
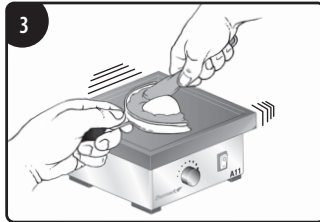
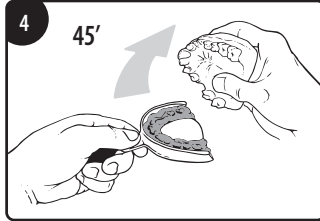
WATER / POWDER RATIO	23 ml / 100 g
WORKING TIME	5'
SETTING TIME (VICAT)*	18'
SETTING EXPANSION 2h*	0.05 %
COMPRESSIVE STRENGTH 1h*	40 MPa

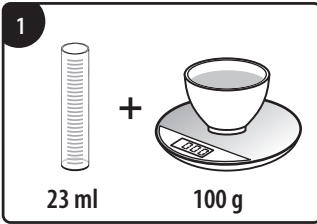
Rapporto acqua/polvere
Water/powder ratio
Wasser/Pulver Verhältnis
Ratio eau/poudre
Dosificación agua/polvo

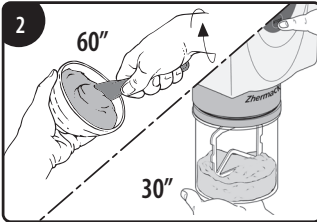Miscelare per 60/30 secondi
Mix for 60/30 seconds
60/30 Sekunden anmischen
Mélanger pendant 60/30 secondes
Mezclar entro 60/30 segundos

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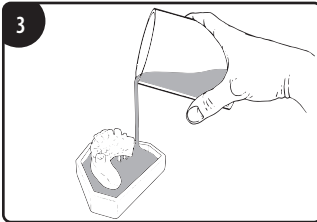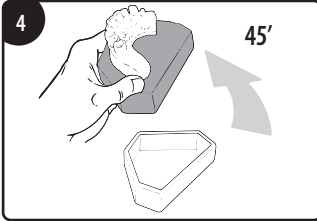EN - WARNING: Be sure to properly disinfect the impression before pouring. PRECAUTIONS: Be aware that during hardening the mixed material may develop enough heat to cause dermal burns. Do not touch the material during hardening. Ensure adequate ventilation. Avoid dust formation and use respiratory and eye protection. Wear gloves to prevent prolonged or repeated skin contact. The times indicated may vary according to room temperature and mixing technique. Dental stone. This product is humidity-sensitive, store in a dry place.

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