



zetalgin

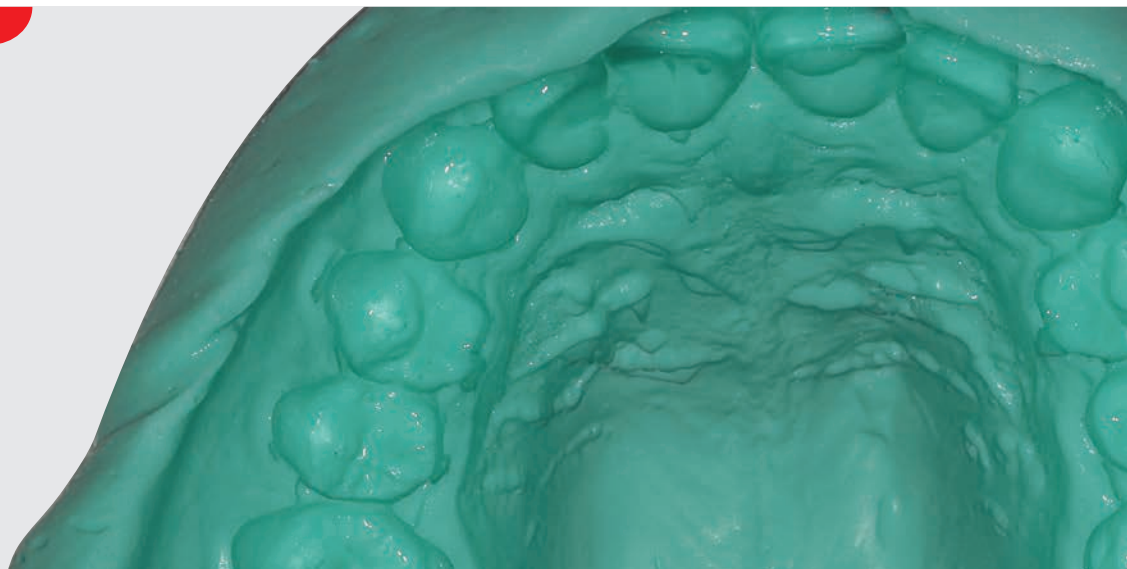

SIMPLE AND EFFECTIVE ALGINATES

Zetalgin alginates are ideal for those professionals searching for simple and effective materials in impression taking:

- ▶ Elasticity for an easy impression removal
- ▶ Easy mixing with fast water absorption for taking impressions quickly
- ▶ 48 hours of dimensional stability for safe model casting after the impression taking
- ▶ Compatibility with mechanical and centrifugal mixers (Zhermack Hurrimix² and Alghamix II) for better efficiency in the dental practice
- ▶ Normal or chromatic version to adapt to the different needs of dental professionals
- ▶ Chromatic version to help professionals by showing the working phase with a change of colour:
 - Purple phase: mixing time
 - Grey phase: working time
 - Green phase: time in mouth

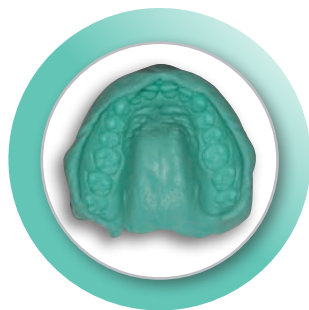

- ▶ A universal preliminary impression alginate to meet the many daily needs of clinical practices



- ▶ Mint scented for both normal and chromatic version



- ▶ Particularly suited to take preliminary impressions to cast both study and antagonist models




- ▶ Chromatic version to assist professionals while taking impressions

Technical details



Product	Time in mouth* (mins)	Strain in compression	Elastic recovery
Zetalgin	1:00	10 %	96 %
Zetalgin chromatic	1:00	10 %	96 %

* the times mentioned must be intended at 23 °C - 73 °F

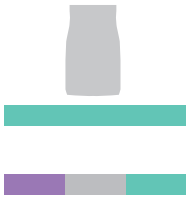

Code	Description
C301001	Zetalgin - 1 bag 453 g
C301004	Zetalgin chromatic - 1 bag 453 g

Accessories

Cod. C300910

Measuring cups set for alginates

Cod. C300992

Rubber mixing bowl

Cod. C300990

Fluorescent coloured alginate spatulas (6 pcs)

Cod. C300960

Long Life Bag, 1 package 100 pcs (bags with hermetic closure to store impressions)